



**A compact
solution to the
complex
problem of
scope storage.**



METRO STARSYS™



METRO STARSYS™

The **STARSYS™** scope storage cabinet is one of the latest Metro® Starsys Modular System solutions for healthcare. Starsys is a fully-integrated, fully-adaptable, modular storage system which includes enclosed storage, procedure carts, transport units, medication carts and WorkCenters. For more information about our scope cabinet, or any of our Starsys solutions, call (800) 433-2232.

Storing expensive, delicate scopes can be a big problem. Fortunately, Metro® has a compact solution – the Starsys™ scope cabinet. With this innovative concept, you can safely store over three scopes per square foot, making it the most storage efficient scope cabinet you can buy.

The Starsys™ scope cabinet is the most convenient cabinet on the market. It gives you easy access to each scope. To better protect your scopes, the cabinet is designed with features such as a solid door, soft vinyl protective coating on the brackets and hoops, and a protective rubber donut on the post.

It all makes for an innovative solution. Which is exactly what people expect from Metro.



**Open
the door
to better
scope
storage.**

Unlike cumbersome tambour style doors, the Starsys scope storage cabinet door opens and closes easily.


Solid door protects scopes. It won't bow inward upon impact.

Vents are strategically located to ensure maximum airflow. Bottom vent is situated near drip pan, which helps prevent odors.



STARSYS

The METRO Carousel™ rack rotates, giving you easy access to any scope.

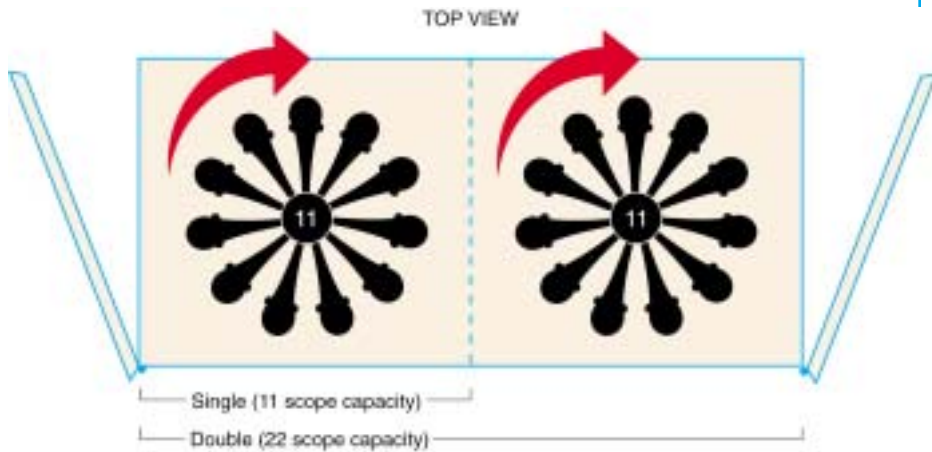

Scope brackets, cord managers and bottom hoop are all adjustable in one-inch (25mm) increments. This allows you to raise or lower them (without tools) to accommodate the scope lengths you are storing.



The mid-level cord manager is used for all types of scopes, and is also orange coated. It keeps cords from tangling and provides strain relief for light boxes or connectors.




Individual scope brackets and cord managers are easy to remove and interchange*. This allows you to mix and match different types of scopes in the same cabinet.



Label holders are available to identify the type of scope or to help locate or track scope inventory by serial numbers.

* Newly developed brackets for changing scope dimensions can be easily retrofitted—saving you the cost of a new cabinet.

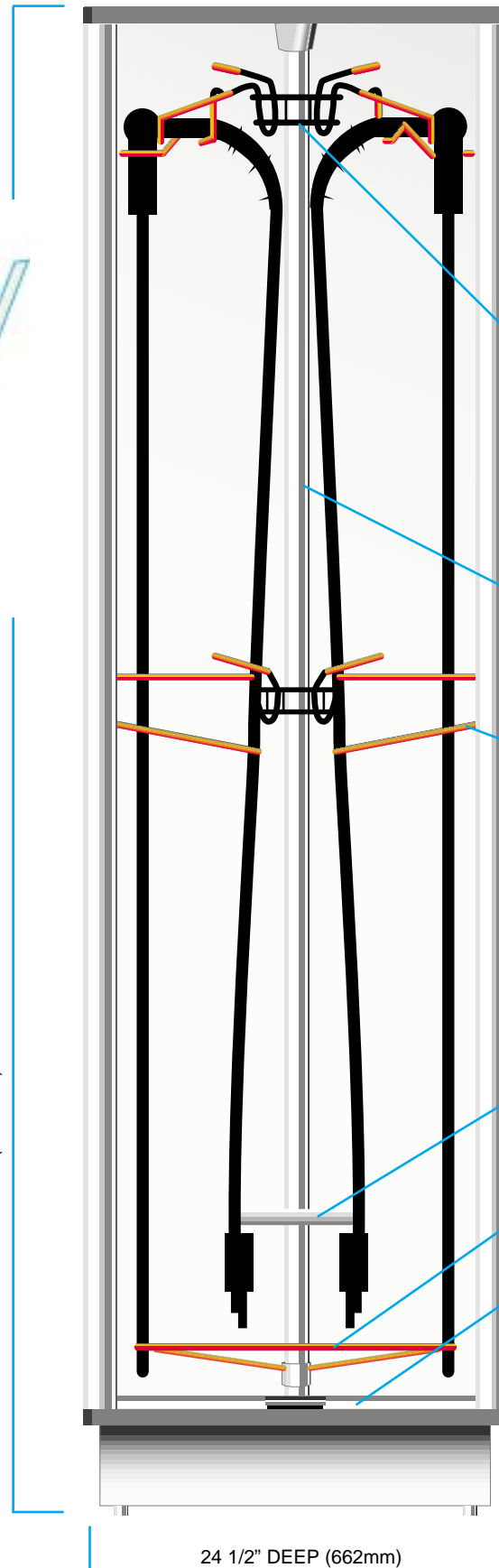

The scope brackets were specifically designed for the different size scopes used. BLUE identifies Pulmonary scopes. ORANGE identifies Gastrointestinal.

Brackets can be easily adjusted up or down depending on the length of scopes being stored.

Carousel post is made of stainless steel to resist corrosion.

Cord manager doubles as turning handle.

Protective rubber donut keeps connector from contacting post.

Soft vinyl coated hoop retains umbilical cords.

Plastic drip pan is easily removed for cleaning and drying.

For more information call 800.433.2232

METRO STARSYS™

www.metro.com



MetroMax® Top Track™ Storage System



Starsys™ Mobile Storage

HUMAN
ENGINEERED
SHELVING,
CARTS AND
WORKCENTERS.



We put space to work.™



Super Adjustable Super Erecta® Shelving System



Starsys™ WorkCenter System

For more information, call 800-433-2232. www.metro.com

